

For patients with osteoarthritis of the knee

50K90 Agilium Softfit

Aside from the indication, the user's day-to-day habits and personal design preferences play a role in selecting the right device for patients with osteoarthritis of the knee. The Agilium Softfit knee orthosis expands the existing portfolio with another useful alternative.

It is a single upright orthosis with a textile base body for patients with unicompartmental osteoarthritis of the knee. Based on the 3-point principle, it is suitable for the treatment of both varus and valgus types. Versions are available for the left and right leg.


The benefits at a glance

- Proven 3-point principle
- Can be fitted to the patient quickly and easily in just a few steps
- Straightforward fitting of varus and valgus types by adjusting the bars
- Good fit without slipping out of place thanks to the innovative material combination



The corrective force of the orthosis can be set with the Allen head screw above the orthosis joint. This makes it easy to adjust the varus/valgus angle of the orthosis.

▶ Further information is available in our guide for patients with osteoarthritis (646D1155=DE)

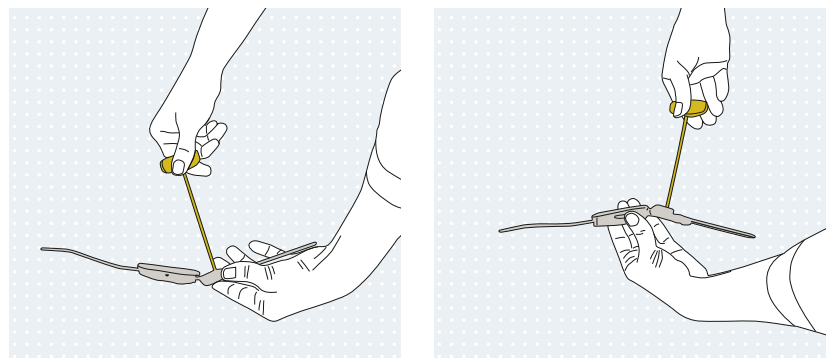

Technical data

		Circumference	S	M	L	XL	XXL
	A	cm	35 – 38	38 – 41	41 – 44	44 – 48	48 – 51
		inches	13.8 – 15.0	15.0 – 16.1	16.1 – 17.3	17.3 – 18.9	18.9 – 20.1
	B	cm	44 – 48	48 – 52	52 – 56	56 – 61	61 – 67
		inches	17.3 – 18.9	18.9 – 20.5	20.5 – 22.0	22.0 – 24.0	24.0 – 26.4